



# ExpressPlex™ 2.0 (384-well) Library Preparation Kit for Illumina® Sequencing Platforms

Catalog numbers:

- 301141:** ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 Reactions - Set A
- 301142:** ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 Reactions - Set B
- 301143:** ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 Reactions - Set C
- 301144:** ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 Reactions - Set D
- 301152:** ExpressPlex 2.0 (384-well) Library Preparation Kit, 1X384 Reactions - Any Index

## User Guide

v20250103

# Table of Contents

<b><u>Introduction</u></b>	3
<b><u>Molecular Diagram</u></b>	4
<b><u>Workflow Diagram</u></b>	5
<b><u>ExpressPlex 2.0 Library Preparation Kit Components</u></b>	6
<b>Considerations before you begin</b>	
<b><u>User-Supplied Reagents, Equipment &amp; Consumables, and Thermal Cycler Programs</u></b>	9
<b><u>Reagent handling</u></b>	
<b>Procedure</b>	
<b><u>ExpressPlex 2.0 Reaction Setup (Using Automated Liquid Handlers)</u></b>	12
<b><u>Library Pooling</u></b>	14
<b><u>Library Pool Purification</u></b>	15
<b><u>Library Quantification and QC</u></b>	17
<b><u>Sequencing on Illumina Platforms</u></b>	20
<b><u>Appendix A: Adjusting the Input Sample Concentration</u></b>	21
<b><u>Appendix B: Sample Sheet and Sequencer Guidelines</u></b>	22
<b><u>Appendix C: ExpressPlex Index Information and Demultiplexing Guidance</u></b>	23
<b><u>Appendix D: Library Quantification Options</u></b>	25

## Introduction




The patent-pending ExpressPlex 2.0 Library Prep Kit comes in both convenient 96-well and 384-well PCR plate configurations for high-throughput multiplexed library preparation. This upgraded version of ExpressPlex uses seqWell's high performance TnX™ transposase that was specifically engineered for NGS library preparation. Amplicons (>350 bp) and plasmid DNA are suitable standard inputs for the kits. ExpressPlex libraries are compatible with the Illumina MiSeq™, NextSeq™, iSeq™, and NovaSeq™ sequencing platforms.

Each ExpressPlex 2.0 (384-well) kit contains sufficient reagents to prepare Illumina-compatible libraries from 384 or 1,536 individual DNA samples. Libraries are typically prepared with 384 samples per pool, up to 1,536 samples per kit. There are four different kits available for preparing libraries from 1,536 samples, providing a total of 6,144 barcode combinations that can be loaded on a single sequencing run.

This multiplexed library preparation procedure is optimized for inputs of 0.5 - 20 ng of plasmid or amplicon DNA per 8 µl reaction, and typically generates library fragment lengths ranging from 400 – 1,200 bp. Library fragment length will depend on the quality and size of input DNA and magnetic bead cleanup ratio used. The primary advantages and benefits of using the ExpressPlex library preparation kits are a streamlined, one-step multiplexed library preparation workflow that automatically normalizes read output per sample across a 40-fold range of DNA input concentration, while minimizing labor and consumable costs. Using the ExpressPlex 2.0 (384-well) kit, a 384-plex library can be prepared for library QC and sequencing in under 120 minutes, with less than 30 minutes of hands-on time.


ExpressPlex 2.0 library preparation kits utilize a proprietary mixture of enzymes to tag input DNA with indexed adapters and amplify libraries all in a single reaction. Each 384-well plate contains 96 different i7-indexed adapters arrayed in a 4x96 grid across the plate. Each grid of 96 i7-indexed adapters contains one i5-indexed adapter, for a total of 384 combinatorial indexed samples per plate. During the single thermal cycling incubation step for 384 samples, the i7-indexed adapters and i5-indexed adapters are simultaneously attached to each of the DNA samples; subsequently, libraries are amplified in segregated amplification reactions. This makes for a highly efficient, one-step multiplexed library prep workflow. Multiple plates of samples may be multiplexed in a single sequencing run, for a total of up to 6,144 samples, provided that the user utilizes a different plate for each set of 384 samples.

## ExpressPlex 2.0 Molecular Diagram

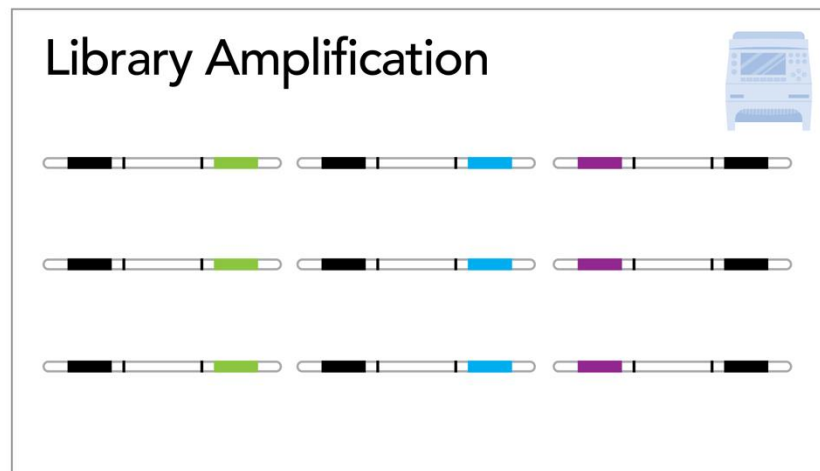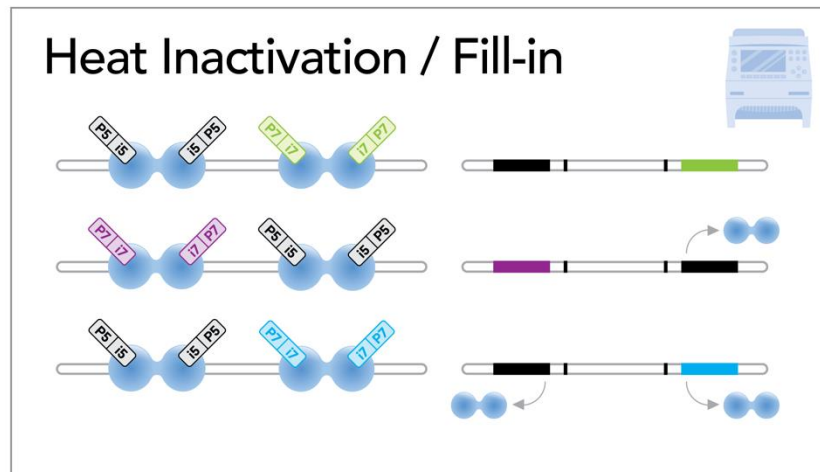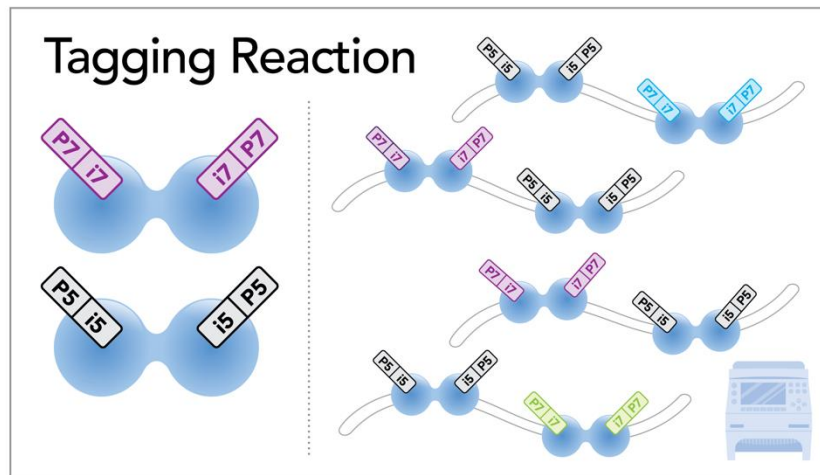

# ExpressPlex 2.0 (384-well)

## Workflow Diagram

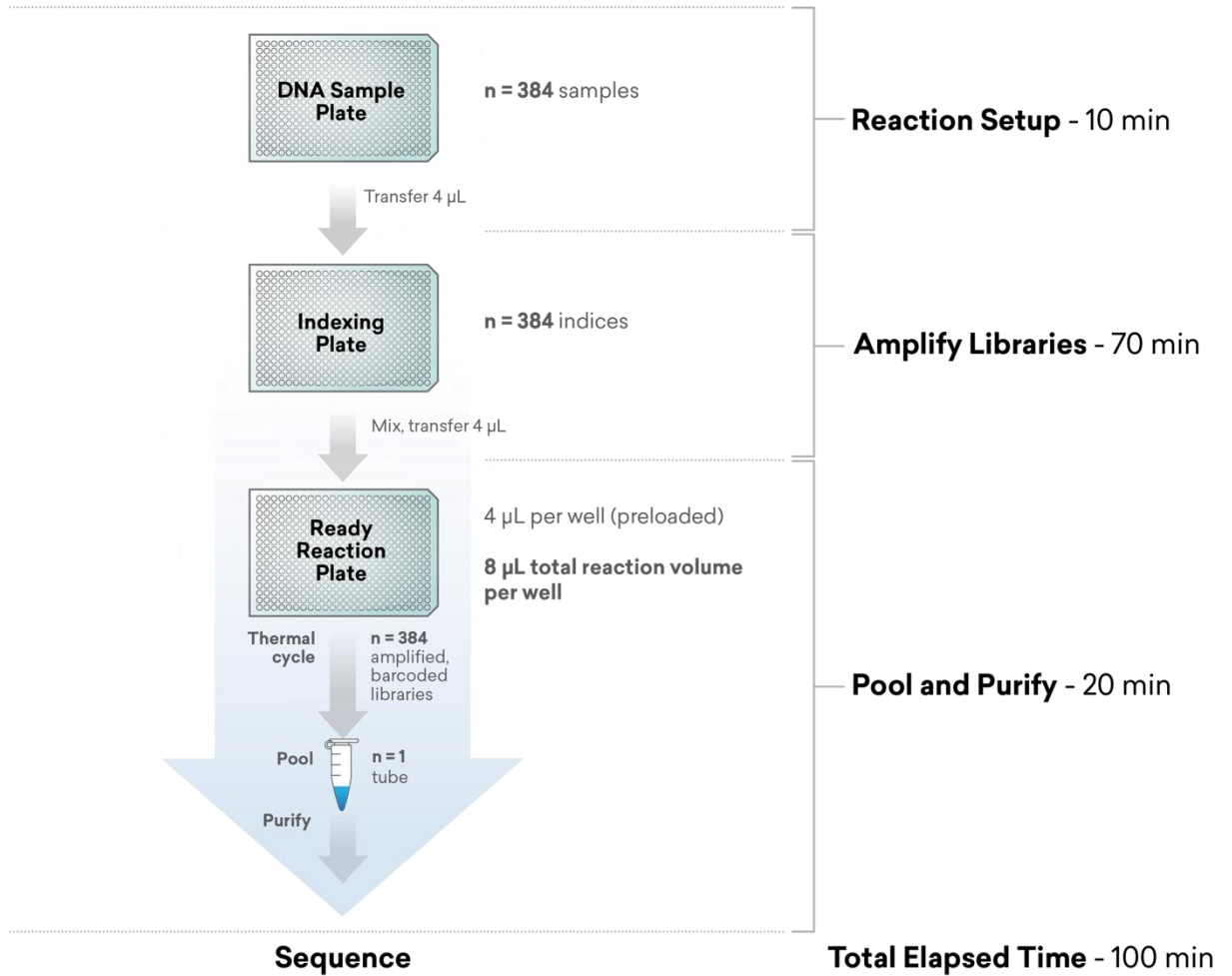

## Kit Components

**Table 1.** ExpressPlex 2.0 (384-well) Library Preparation Kit Components (4X384 Reactions)

### ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 - Set A

Catalog No.: 301141

Item	Component	REF	Description	Storage	Qty
1	Box 1 - 4X384 Indexing Reaction Plates – A-1000, A-2000, A-3000, A-1000	301147	Indexing Reagent Plates (384-well) in a full-skirt, colored PCR plates	-20° C	4
2	Box 2 - 4X384 Ready Reaction Mix Plates	301151	Ready Reaction Plates (384-well) in a full-skirt, clear PCR plates	-20° C	4
3	MAGwise paramagnetic beads*	101003	5 ml of MAGwise beads in a 10 ml screwcap tube	4° C	1

### ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 - Set B

Catalog No.: 301142

Item	Component	REF	Description	Storage	Qty
1	Box 1 - 4X384 Indexing Reaction Plates – B-1000, B-2000, B-3000, B-4000	301148	Indexing Reagent Plates (384-well) in a full-skirt, colored PCR plates	-20° C	4
2	Box 2 - 4X384 Ready Reaction Mix Plates	301151	Ready Reaction Plates (384-well) in a full-skirt, clear PCR plates	-20° C	4
3	MAGwise paramagnetic beads*	101003	5 ml of MAGwise beads in a 10 ml screwcap tube	4° C	1

**ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 - Set C****Catalog No.: 301143**

Item	Component	REF	Description	Storage	Qty
1	Box 1 - 4X384 Indexing Reaction Plates – C-1000, C-2000, C-3000, C-4000	301149	Indexing Reagent Plates (384-well) in a full-skirt, colored PCR plates	-20° C	4
2	Box 2 - 4X384 Ready Reaction Mix Plates	301151	Ready Reaction Plates (384-well) in a full-skirt, clear PCR plates	-20° C	4
3	MAGwise paramagnetic beads*	101003	5 ml of MAGwise beads in a 10 ml screwcap tube	4° C	1

**ExpressPlex 2.0 (384-well) Library Preparation Kit, 4X384 - Set D****Catalog No.: 301144**

Item	Component	REF	Description	Storage	Qty
1	Box 1 - 4X384 Indexing Reaction Plates – D-1000, D-2000, D-3000, D-4000	301150	Indexing Reagent Plates (384-well) in a full-skirt, colored PCR plates	-20° C	4
2	Box 2 - 4X384 Ready Reaction Mix Plates	301151	Ready Reaction Plates (384-well) in a full-skirt, clear PCR plates	-20° C	4
3	MAGwise paramagnetic beads*	101003	5 ml of MAGwise beads in a 10 ml screwcap tube	4° C	1

\*A larger 15ml bottle of MAGWise paramagnetic beads is available (REF 101002) for purchase if desired for easier automation.

**Table 2.** ExpressPlex 2.0 (384-well) Library Preparation Kit Components (1X384 Reactions)**ExpressPlex 2.0 (384-well) Library Preparation Kit, 1X384 – Any Index****Catalog No. 301152**

<b>Item</b>	<b>Component</b>	<b>REF</b>	<b>Description</b>	<b>Storage</b>	<b>Qty</b>
1	Indexing Reaction Plate (any index); Ready Reaction Mix Plate	301153	Indexing Reagent Plate (384-well) in a full-skirt, colored PCR plate; Ready Reaction Plate (384-well) in a full-skirt, clear PCR plate	-20° C	1
2	MAGwise paramagnetic beads*	101003	5 ml of MAGwise beads in a 10 ml screwcap tube	4° C	1



# User-Supplied Reagents, Equipment, Consumables, and Thermal Cycler Program

## **Reagents**

- 80% Ethanol (molecular biology grade)
- Tris-HCl, pH 8.0
- Ultrapure Water (PCR grade)
- Quant-iT™ PicoGreen™ dsDNA Assay Kits (ThermoFisher P/N: P7589), or other validated dsDNA quantification assay
- Electrophoretic analysis reagents - Agilent TapeStation (High Sensitivity D5000 kit), Bioanalyzer (High Sensitivity DNA or DNA7500 kits), or Fragment Analyzer (High Sensitivity NGS Fragment Analysis Kit, DNF-474)

## **Equipment & Consumables**

- Single-channel pipettors (1-20 µl, 20-200 µl, 100-1,000 µl)
- Multichannel pipettors (1-10 µl, 10-200 µl)
- Conventional liquid handler with a 96- or 384-tip pipetting head
- Pipette tips (low-retention barrier tips)
- Eppendorf Tubes® (1.5 ml and 2 ml, DNA LoBind Tubes)
- PCR plate seals (must be evaporation-resistant)
- 384-well thermal cycler (compatible with fully skirted Bio-Rad Hard-Shell® 384-Well PCR Plates )
- Magnetic stand for 1.5 ml tubes
- 0.2 ml 8-tube PCR strips and caps/seals
- Benchtop centrifuge to pulse-spin tubes and 8-tube PCR strips
- Plate centrifuge
- Vortex mixer
- Electrophoretic analysis equipment - Agilent TapeStation, Bioanalyzer, or Fragment Analyzer
- Fluorometer for dsDNA quantification assay

**Thermal Cycler Program EP2\_TAG\_AMP (all with lid-heating on at 105°C):**

15	min	55°C	Tagging Reaction
5	min	75°C	Heat Inactivation / Fill-in
3	min	68°C	
5	min	79°C	
3	min	68°C	
5	min	83°C	
1	min	98°C	Initial denaturation
15	sec	98°C	Library Amplification (PCR) <b>Use 12 cycles for plasmids</b> <b>Use 15 cycles for amplicons</b>
30	sec	64°C	
1	min	72°C	
5	min	72°C	Final extension
Hold		4°C	

## Before starting the procedure:

**Measure and adjust input DNA concentration.** It is recommended that the user assay the DNA concentration of each plate of samples to be processed by PicoGreen or other validated dsDNA assay. This can be accomplished by measuring the concentration of a representative selection of samples. Globally adjust the average concentration of input DNA across each plate to 1.25 ng/ $\mu$ l (2.5 ng calculated input into final reaction) in 10 mM Tris-HCl, pH 8.0. Do not use EDTA-containing solutions (e.g., TE buffer) to dissolve or dilute input DNA because EDTA can suppress enzymatic activity. *Refer to Appendix A for more detailed information on globally adjusting the average input DNA concentration.*

**Program liquid handler.** For assistance, contact support@seqwell.com

**Program thermal cycler.** For convenience, set up the thermal cycler program described in the previous page before starting.

**Pulse-spin kit components.** Liquids can condense or shift locations inside containers during shipment or storage. Before using the reagent plates and tubes, pulse-spin in a suitable centrifuge to gather the reagents at the bottom of the well or tube. If the kit components freeze during shipment or storage, thaw, mix and pulse-spin before use.

**Equilibrate MAGwise Paramagnetic Beads to room temperature.** MAGwise beads can be stored for up to 2 weeks at room temperature or for longer periods at 2 - 8°C. If stored cold, warm at room temperature for 30 minutes before use. Vortex to thoroughly resuspend beads prior to use. To transfer volumes accurately, pipette slowly and do not pre-wet pipette tips.

**Prepare 80% ethanol fresh daily.** You will need ~5 ml per sample plate.

**Prepare 10 mM Tris-HCl, pH 8.0.** Prepare 10 mM Tris-HCl, pH 8.0 from a concentrated stock solution diluted with ultrapure water (both molecular-biology grade). Do not use EDTA-containing solutions (e.g., TE). User will need ~50  $\mu$ l per 384-well sample plate.

**Safe stopping points** are indicated in the protocol. For optimal results, proceed directly to the next step unless a safe stopping point is indicated.

# Procedure

## 1. ExpressPlex Reaction Setup

### **Note: Key Items Before Starting:**

- Performing setup via automated liquid handler is highly recommended.
- Refer to **Appendix A** for instructions to dilute input samples.
- If processing more than one plate, multiple reaction plates may be set up and thermal cycled at the same time before proceeding to subsequent steps. (Immediate continuation after cycling is recommended)
- User will pipette 4  $\mu$ l of DNA into the tagging reagent, resulting in a 1:1 dilution of the DNA sample prior to beginning the next step.

ExpressPlex reactions can be set up at room temperature. Performing setup via automated liquid handler is highly recommended.

- a. Pulse-spin one **Indexing Reagent Plate** and one **Ready Reaction Plate** in a tabletop centrifuge for 1 minute at 1000 x g.
- b. Pre-label each **Ready Reaction Plate** with the number from the **Indexing Reagent Plate** and a name to allow easy identification of your samples.
- c. Transfer 4  $\mu$ l of input DNA to each well (one sample per well) of the **Indexing Reagent Plate**. Mix thoroughly by pipetting up and down ten times at 4  $\mu$ l, being careful not to introduce excessive bubbles.
- d. Using the same tips, immediately transfer 4  $\mu$ l of the resulting **DNA/Indexing Reagent mixture** to each corresponding well of the **Ready Reaction Plate**. Mix thoroughly by pipetting up and down ten times at 4  $\mu$ l, being careful not to introduce excessive bubbles. Note: 4  $\mu$ L of unused volume will remain behind in the Indexing Plate after transfer. This volume should be discarded.

*Optional:* If excessive bubbles are present after reaction set-up, use a tabletop centrifuge to collapse the bubbles prior to proceeding.

- e. Seal the reaction plate, pulse-spin, then transfer to a thermal cycler, and run the **EP2\_TAG\_AMP** program below, with lid-heating on at 105°C:

15	min	55°C	Tagging Reaction
5	min	75°C	Heat Inactivation / Fill-in
3	min	68°C	
5	min	79°C	
3	min	68°C	
5	min	83°C	
1	min	98°C	Initial denaturation
15	sec	98°C	Library Amplification (PCR) <b>Use 12 cycles for plasmids</b> <b>Use 15 cycles for amplicons</b>
30	sec	64°C	
1	min	72°C	
5	min	72°C	Final extension
Hold		4°C	

***SAFE STOPPING POINT***

**Proceed immediately to the next step or  
store the amplified libraries at -20°C.**

## 2. Library pooling

- a. After library amplification, pulse-spin the **Ready Reaction Plate**, and then remove the plate seal.
- b. Using a liquid handler, pool 4  $\mu\text{l}$  of each amplified library from each 96-well quadrant into a single pre-labeled 96-well PCR plate. The four quadrants will be combined, 4 samples per well, in this plate.

*Note:* The same pipette tips may be reused for pooling because the samples are fully barcoded and amplified.

- c. Using 10  $\mu\text{l}$  from each combined well, proceed to pooling all samples, by column, into a separate plate or 8-tube PCR strip.
- d. After mixing by pipetting, transfer 100  $\mu\text{l}$  from each well of the plate or 8-tube PCR strip(s) into a pre-labeled 2 ml LoBind tube. The total volume should be approximately 800  $\mu\text{l}$ . Approximately 20  $\mu\text{l}$  will be left behind in each well or strip tube and may be discarded.

*Optional:* If excessive bubbles are present after pooling, use a tabletop centrifuge to collapse the bubbles at 1,000 x g for 1 minute prior to proceeding. Also, you may freeze any unpurified amplified libraries remaining in the Ready Reaction Plate, providing an option to purify more library later if any sample(s) should require additional sequencing depth.

### 3. Library Pool Purification

- a. Vortex (or vigorously pipet) room temperature **MAGwise** to ensure that the beads are fully resuspended before use.
- b. Add 0.75 volumetric equivalent of **MAGwise** to the pooled ExpressPlex library (e.g., add 600  $\mu$ l of MAGwise to 800  $\mu$ l of pooled library), and mix thoroughly by pipetting up and down  $\geq 10X$ .

*Note:* Use 1 volumetric equivalent of **MAGwise** for libraries derived from amplicons in order to recover the shorter fragments.

- c. Incubate on the bench for 5 minutes to allow the DNA to bind.
- d. Place the tube on a magnetic stand and let the bead pellet form on the inner wall of the tube and wait until the supernatant appears completely clear (approximately 5 minutes).
- e. Remove and discard the supernatant with pipette. Be careful not to disturb the bead pellet.
- f. Wash beads with 80% ethanol.
  - i. With the tube still in the magnetic stand, add 1.5 mL 80% ethanol to cover the bead pellet.
  - ii. After  $\geq 30$  seconds, remove and discard the supernatant being careful to not dislodge the bead pellet.
- g. Repeat the previous step for a total of 2 washes with 80% ethanol. [**Do not air dry the bead pellets – proceed immediately to the next step**]

*Useful tip:* After using a large pipette tip to remove the waste ethanol from the second wash, briefly pulse-spin, and then use a smaller pipette tip to remove any residual volume, if visible.

- h. Immediately remove the tube from the magnetic stand and pipette 32  $\mu$ l of 10 mM Tris-HCl, pH 8 on top of the bead pellet. Pipette the solution along the inner wall of the tube multiple times to thoroughly resuspend the bead pellet.
- i. Incubate the tube on the bench for at least 5 minutes to elute the libraries from the beads.
- j. Return tube to magnetic stand and allow the bead pellet to reform on the inner wall of the tube (~2 minutes).

- k. When the supernatant has cleared completely, carefully transfer 30  $\mu$ l of DNA eluate to a clean 1.5 ml LoBind tube. The transferred eluate contains the purified pooled library.

***SAFE STOPPING POINT***

**Proceed immediately to the next step  
or store the pooled purified library at -20°C.**




## Library Quantification and QC


**Library quantification and QC with electrophoretic analysis:** Run the pooled, purified ExpressPlex 2.0 library on the Agilent TapeStation (High Sensitivity D5000), Bioanalyzer (High Sensitivity DNA or DNA7500 kits), or Fragment Analyzer (High Sensitivity NGS Fragment Analysis Kit, DNF-474) following the manufacturer's instructions for these instruments. For optimal sequencing results with ExpressPlex 2.0 plasmid libraries, use a region analysis for fragments of **200 - 1,200 bp** to determine the average clusterable fragment length for size adjustment. For amplicon libraries, adjust the region size based on the input amplicon length to exclude unfragmented amplicons.

If using electrophoretic analysis for library quantification, use the region analysis function and calculate the library concentration from the clusterable fragment region. Alternatively, use Quant-iT PicoGreen dsDNA Assay, Qubit dsDNA HS Assay (see below) or other validated dsDNA quantification assay to quantify the library pool, and calculate amount of clusterable fragments by multiplying by the percentage of sample in the region analysis (200-1200 bp) for sequencer loading. For a more accurate measurement of the concentration of clusterable library fragments, use qPCR (see below).

Refer to the figures below for representative trace for pooled, purified libraries run on the TapeStation. Refer to *Appendix D* for quantification methods by Agilent TapeStation for ExpressPlex 2.0 libraries.

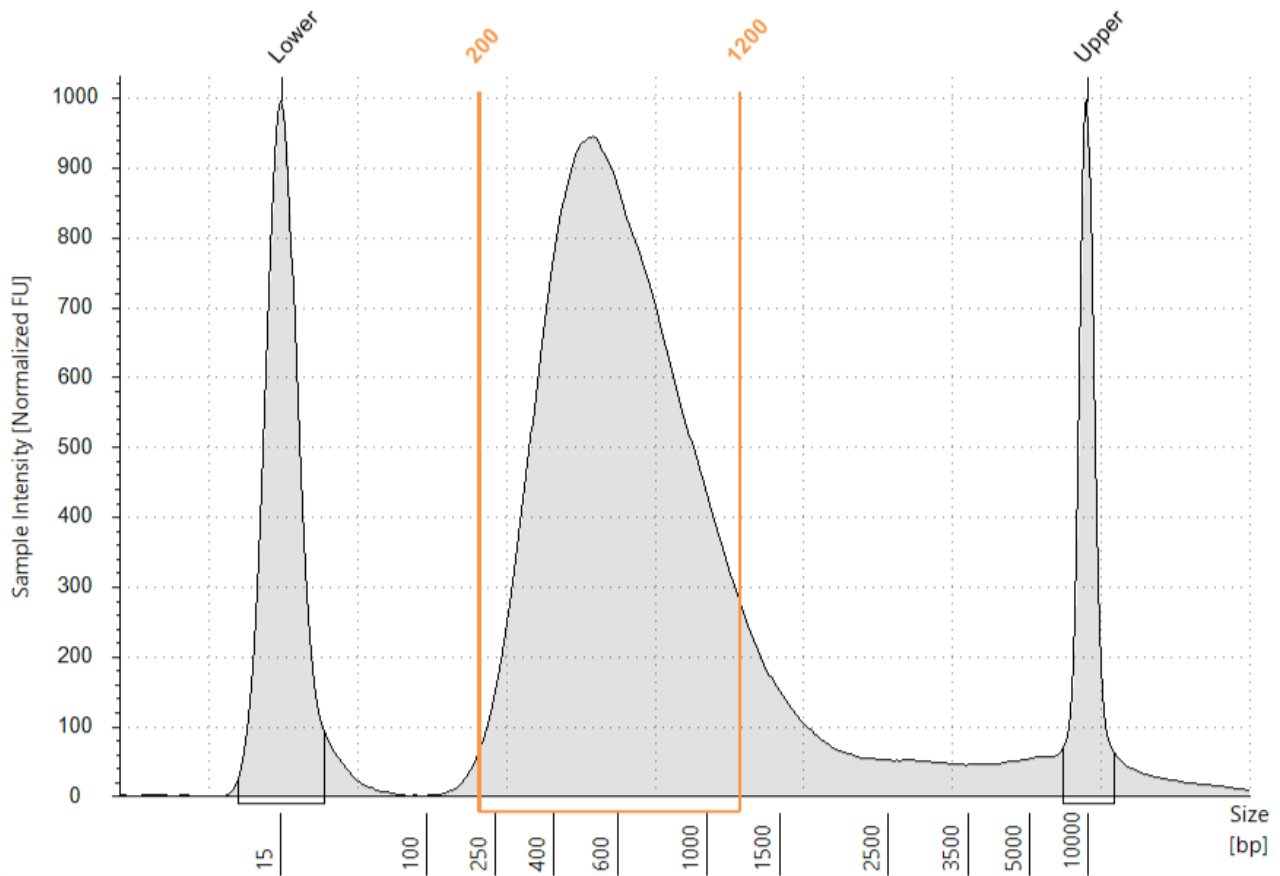

**Figure 1.** Representative TapeStation electropherogram of a pooled, purified ExpressPlex 2.0 plasmid library (diluted 20-fold prior to electrophoresis). The analysis region (orange) shows the range of clusterable fragments.

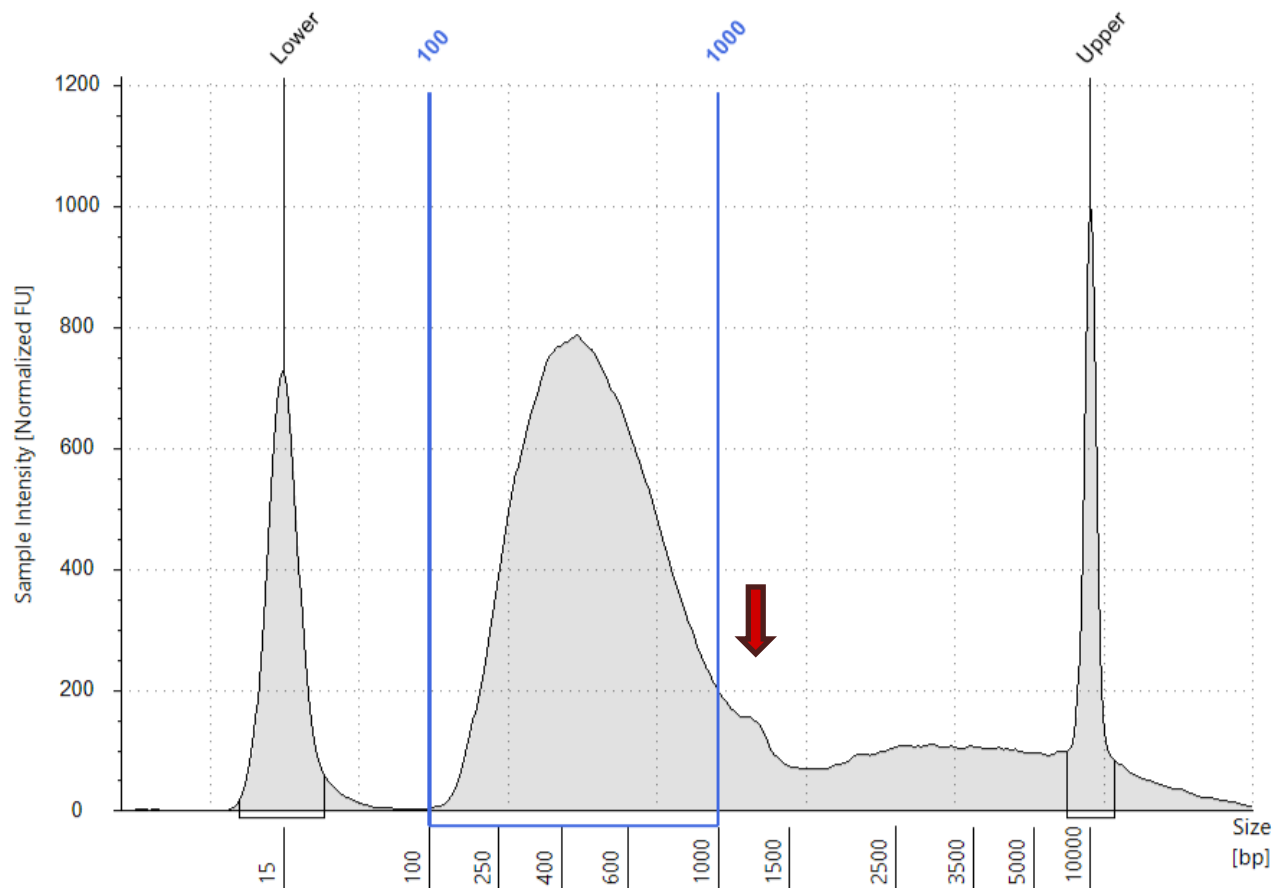

**Figure 2.** Representative TapeStation electropherogram of a pooled, purified ExpressPlex 2.0 amplicon library (diluted 20-fold prior to electrophoresis). The analysis region (blue) shows the range of clusterable fragments. Arrow (red) shows a small amount of unfragmented amplicon DNA.

**Library quantification with dsDNA specific fluorometric method:** ExpressPlex 2.0 libraries are double stranded. Use 2 µL to quantify the pooled, purified library using the Qubit dsDNA assay, Quant-iT PicoGreen dsDNA assay or other validated dsDNA quantification assay. Follow the manufacturers' instructions for the appropriate dilutions. Use the average clusterable fragment size as determined by electrophoresis to calculate the library concentration and the dilution factor required for loading your sequencing system.

**Library quantification with qPCR assay:** qPCR is a very sensitive method of quantifying library fragments that have both adaptor sequences on either end which will subsequently form clusters on a flow cell. Use 2 µl of the purified, multiplexed ExpressPlex 2.0 library for qPCR analysis with Illumina qPCR primer 1.1/2.1 described on Illumina's website here ([https://support.illumina.com/content/dam/illumina-support/documents/documentation/chemistry\\_documentation/qpcr/sequencing-library-qpcr-quantification-guide-11322363-c.pdf](https://support.illumina.com/content/dam/illumina-support/documents/documentation/chemistry_documentation/qpcr/sequencing-library-qpcr-quantification-guide-11322363-c.pdf)). Follow kit and instrument instructions for appropriate conditions and dilutions.

## Sequencing on Illumina platforms

**Read configuration:** ExpressPlex Library Prep kit libraries are dual indexed with 10 base indexes. Index lists can be found in Appendix C and in spreadsheet form on the seqWell [website](#). ExpressPlex libraries are sequenced using the same primers as standard Illumina libraries, and consequently, custom sequencing primers are not needed. Longer reads deliver greater read depth and read lengths of 2 x 150 or greater are recommended for *de novo* assembly. The i7 and i5 index reads must be 10 bases long for ExpressPlex libraries, although the index reads and non-index reads can be adjusted for different sequencing kits, speed, or read depth requirements.

**Library dilution, denature and sequencing:** Follow Illumina's guidelines for appropriate dilution to the starting concentration and loading procedures specific to your sequencing system. Not all DNA fragments can efficiently generate clusters. Longer library fragments generally do not cluster as efficiently on a cell. It is essential to optimize the loading concentration based on your sample type, library QC methods and library fragment distribution. Try adjusting loading concentration over subsequent sequencing runs to optimize cluster density and sequencing run performance.

## Appendix A: Adjusting the Input Sample Concentration

ExpressPlex 2.0 (384-well) library prep kits perform optimally with 0.5 - 20 ng of total dsDNA input per reaction, and individual adjustment of each sample concentration is not necessary.

If the method used to purify DNA upstream of library prep is well-characterized and generates relatively consistent concentrations of DNA per sample, it may be sufficient to only assay several samples from the DNA sample plate (*i.e.*, spot-check the dsDNA concentration using a PicoGreen Assay to estimate the average DNA concentration across all samples).

If all of the samples already fall within the 0.25 – 10 ng/μl range, no adjustment is required. If, however, the average concentration of all of the samples exceeds 5 ng/μl, calculate the global dilution factor using the formula below:

$$\text{Global dilution factor (X)} = \frac{\text{Average assayed dsDNA concentration (ng/}\mu\text{l)}}{1.25 \text{ ng/}\mu\text{l}}$$

The global dilution factor is applied to the input samples so that the average DNA concentration across all samples will be approximately 1.25 ng/μl (*i.e.*, resulting in an average of 5 ng of input being added per full ExpressPlex reaction volume).

**Note:** If the DNA concentration of the input samples is not easily confined to an approximately 40-fold range (maximum concentration divided by minimum concentration), or, if an average sample concentration of  $\geq 0.25$  ng/μl cannot be routinely achieved, consider optimizing the method used to generate input DNA.

**Important Reminder:** Do not use EDTA-containing solutions (*e.g.*, TE buffer) to dissolve or dilute input DNA because EDTA can suppress enzymatic activity.

## Appendix B: Sample Sheet and Sequencer Guidelines

Illumina sequencing systems and chemistries differ in their use of sample sheets, availability of on-instrument demultiplexing, optimal color balancing, and run setup methods. If you have questions for your specific sequencer, contact Illumina tech support for the best guidance on setting up your run using your sequencing platform.

ExpressPlex 2.0 libraries are similar to the Nextera Adapter sequences and are dual-indexed using 10 base indexes for both the i7 and i5 index sequences. ExpressPlex libraries do NOT require custom sequencing primers.

All Illumina sequencers read the i7 index in the forward direction (as listed in Appendix C). The i5 index, however, is read differently on different sequencers depending on the version/chemistry of the sequencing kits. In this case, if using an Illumina sample sheet template, enter the i5 index in the forward direction as the sample sheet will auto-generate the reverse complement if needed. If demultiplexing using bcl2fastq and using a sequencer that reads the reverse complement of the i5 index sequence, enter the reverse complement of the i5 index shown in the ExpressPlex 2.0 index list on seqWell's website: <https://seqwell.com/expressplex-2-0-library-prep-kit/>

Please refer to Illumina's website for most up to date index sequencing guidelines:

<https://support.illumina.com/downloads/indexed-sequencing-overview-15057455.html>

## Appendix C: ExpressPlex Index Information and Demultiplexing Guidance

ExpressPlex utilizes a combinatorial dual indexing (CDI) strategy. Each 384-well plate contains 96 different i7-indexed adapters, each placed in 4 adjacent wells. Each placement of an individual i7-indexed adapter is combined with 4 different i5-indexed adapters, resulting in 384 CDI combinations per plate. See figure 4 below. All plates from set A (plates A-1000, A-2000, A-3000, and A-4000) contain the same 96 i7's, but utilize 16 different i5's (4 per plate) across the set to achieve 1,536 different index combinations. Each subsequent set (B, C, or D) contains a different set of 96 i7-indexed adapters than any other set but repeats the 16 different i5-indexed adapters. If using all four sets (16 plates), 6,144 total index combinations are available to multiplex into a single run. Please refer to the ExpressPlex index list in spreadsheet format on seqWell's website: <https://seqwell.com/expressplex-2-0-library-prep-kit/>, for a complete list of all i7 and i5 indexes.

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24
A	1	1	9	9	17	17	25	25	33	33	41	41	49	49	57	57	65	65	73	73	81	81	89	89
B	1	1	9	9	17	17	25	25	33	33	41	41	49	49	57	57	65	65	73	73	81	81	89	89
C	2	2	10	10	18	18	26	26	34	34	42	42	50	50	58	58	66	66	74	74	82	82	90	90
D	2	2	10	10	18	18	26	26	34	34	42	42	50	50	58	58	66	66	74	74	82	82	90	90
E	3	3	11	11	19	19	27	27	35	35	43	43	51	51	59	59	67	67	75	75	83	83	91	91
F	3	3	11	11	19	19	27	27	35	35	43	43	51	51	59	59	67	67	75	75	83	83	91	91
G	4	4	12	12	20	20	28	28	36	36	44	44	52	52	60	60	68	68	76	76	84	84	92	92
H	4	4	12	12	20	20	28	28	36	36	44	44	52	52	60	60	68	68	76	76	84	84	92	92
I	5	5	13	13	21	21	29	29	37	37	45	45	53	53	61	61	69	69	77	77	85	85	93	93
J	5	5	13	13	21	21	29	29	37	37	45	45	53	53	61	61	69	69	77	77	85	85	93	93
K	6	6	14	14	22	22	30	30	38	38	46	46	54	54	62	62	70	70	78	78	86	86	94	94
L	6	6	14	14	22	22	30	30	38	38	46	46	54	54	62	62	70	70	78	78	86	86	94	94
M	7	7	15	15	23	23	31	31	39	39	47	47	55	55	63	63	71	71	79	79	87	87	95	95
N	7	7	15	15	23	23	31	31	39	39	47	47	55	55	63	63	71	71	79	79	87	87	95	95
O	8	8	16	16	24	24	32	32	40	40	48	48	56	56	64	64	72	72	80	80	88	88	96	96
P	8	8	16	16	24	24	32	32	40	40	48	48	56	56	64	64	72	72	80	80	88	88	96	96

Numbers 1-96 = i7 location

i5-1	i5-2
i5-3	i5-4

**Figure 3.** Representative layout of the i7- and i5-indexed adapters on an individual plate.

**Optimal Demultiplexing Guidance:**

In accordance with Illumina's guidelines, it is recommended to allow for 1 mismatch in the barcode reads during demultiplexing. This allows for the recovery of substantially more reads than 0 mismatch at a reduced risk, due to sufficient Hamming distances between barcodes.

Refer to *Illumina's mismatch guidelines here*: [https://knowledge.illumina.com/software/general/software-general-reference\\_material-list/000007484](https://knowledge.illumina.com/software/general/software-general-reference_material-list/000007484)

For the XLEAP chemistries:


For NovaSeq X, NextSeq 2000, and MiSeqi100 running XLEAP-SBS chemistry, an increased rate of 1 mismatch in the barcode reads has been observed. Hence, it is especially important to allow for 1 mismatch on these instruments to recover all expected data for each barcode.



## Appendix D: Library Quantification Options

In this example, Agilent TapeStation was employed to quantify a pooled, purified plasmid ExpressPlex 2.0 library for sequencer loading.

Refer to the figure below for library QC results.

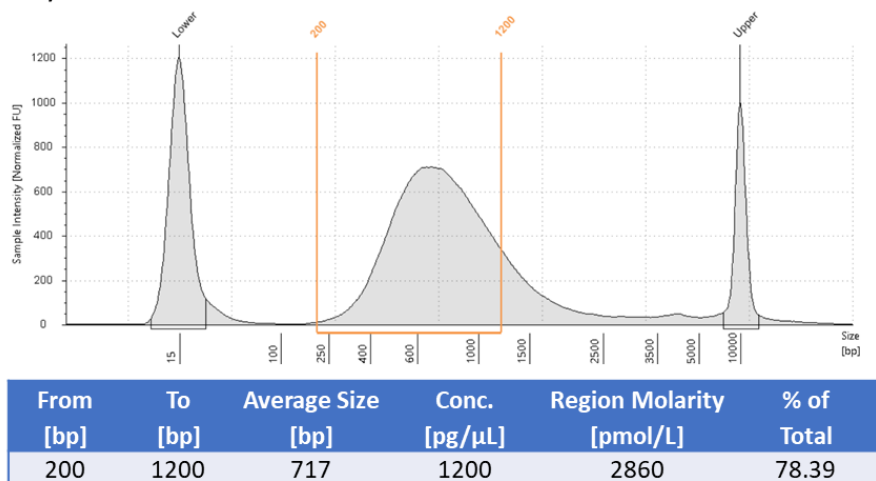

**Figure 4.** Example of library quantification results for a pooled, purified plasmid library by Library QC on Agilent TapeStation.

Use the following formula to convert library concentration from ng/μl to nM:

$$\text{Library concentration [nM]} = \frac{\text{Library concentration [ng/μL]}}{660[\text{g/mol}] \times \text{Average fragment size [bp]}} \times 10^6$$

Library quantification by TapeStation **ONLY** (Library diluted 20-fold prior to electrophoresis):

- Clusterable region average library size: 717 bp
- Clusterable region library concentration: 1200 pg/μl × 20 = 24.0 ng/μl
- Estimated library concentration:  $\frac{24.0 [\text{ng/μl}]}{660[\text{g/mol}] \times 717 [\text{bp}]} \times 10^6 = 50.7 \text{ nM}$

Version	Release Date	Prior Version	Description of changes
V20250103	January 3, 2025	N/A	First version

## Technical Assistance

For technical assistance with ExpressPlex 2.0 Library Preparation Kit, contact seqWell Technical Support.

E-mail: [support@seqwell.com](mailto:support@seqwell.com)

seqWell, Inc.  
66 Cherry Hill Drive  
Beverly, MA 01915  
USA  
+1-855-737-9355  
<https://seqwell.com/>

© 2025 seqWell. All rights reserved.  
For Research Use Only – Not for Diagnostic Use  
UserGuide\_EP38420250103

seqWell, the seqWell logo, TnX, and ExpressPlex are trademarks of seqWell, Inc. All other trademarks are property of their respective owners.